IMPLANTOLOGIE

IMPLANTOLOGY
IMPLANTOLOGIE
IMPLANTOLOGÌA

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Implantological hammer / Marteau pour l'implantologie / Martillo para implantologìa

1795
25cm
Hammer für Knochenschlagdose
OMBREDANNE
rostfrei, Gewicht: 720g
Hammer acc. to Ombredanne
for bone centrifuge, rustfree, weight: 720g
Maillet d'après Ombredanne pour boîte de broyage des os inoxydable, poids: 720g
Martillo según Ombredanne para recipiente de martilleo de hueso inoxidable. peso: 720g

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Anatomical tweezers with tungsten carbide (TC) inserts for retracting / gripping thin tissue flaps atraumatically or grasping membranes firmly while cutting

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This instrument allows rapid and measured bone removal close to the insertion site. The high-quality hardened blades enable bone to be removed gently and deliberately in thin layers without excessive force. To maintain optimal sharpness, we recommend grinding the blades with an Arkansas stone (CM no. 995) after each use (prior to sterilisation).

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
BONE SCRAPER
- fast and simple generation of autologous bone materials
- the blade enables the generation of cortical and spongeous bones
- the recovered bone is collected in the chamber and mixed with blood
- one can in the norm do without tool aids such as bone filter, drilling awl or saws
- ergonomically perfect design

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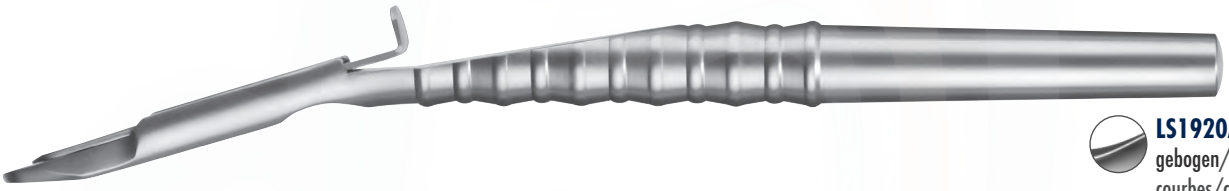Instrument for introducing bone material
- The bone can be removed from the bone mixer tin quickly and safely with the curette. It is then applied to the appropriate site.

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
Bone Mill, rustfree with Titanium millworks
- the Titanium jaws are exchangable
- high stability and easy handling
- effortless cleaning

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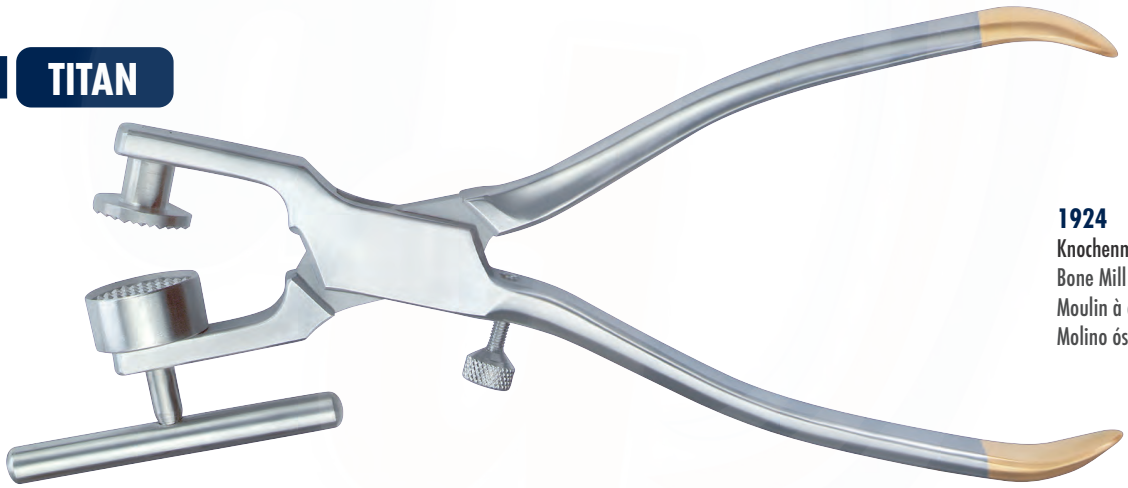
Bone Percussion Tin
- for longer pieces of bone
- all Titanium to insure an optimal tissue compatibility (biocompatible).

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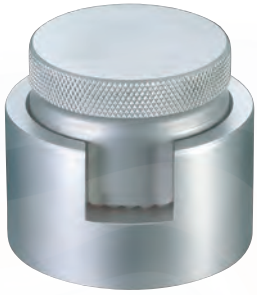- Sinus-lifting: Describes an augmentation method in the immediate vicinity beneath the sinus in cases with atrophied mandibles. Sinus-lifting involves first releasing the mucous membrane of the nose from the bone and filling the remaining space with autologous bone
- Instruments 1142/1 and 2 are for raising and simultaneously displacing the sinus membrane cranio-medially in the sinus

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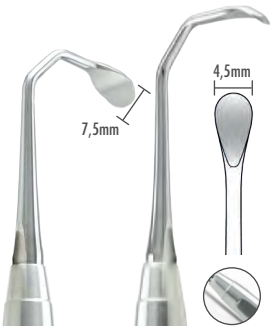
Marking as safety zone
available in different widths

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Marking as safety zone available in different widths, conical, convex working end

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Marking as safety zone available in different widths cylindrical, concave and medium sharp working tip

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- straight or curved design (straight = maxilla, curved = mandible)
- different widths from 3 to 8mm, recommendation for single tooth implants: up to 4mm width, larger jaw segments: up to 8mm width
- depth markings on the osteotome (black region = safety marking above 10mm)
- recommendation for use: min. ridge width of 2-3mm, min. bone height on 10mm, max. splitting depth up to 70 percent of total bone height

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
- straight or bayonet design (straight = for anterior and posterior tooth regions) (curved = for molar region)
- different widths with cylindrical working end
- the working tip is concave and has a medium sharp edge, which allows the bone to be pushed under the sinus mucosa and compressed

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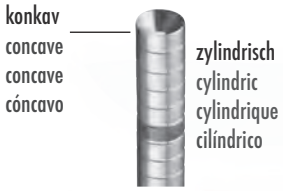
- straight or bayonet design (straight = for anterior and posterior tooth regions) (curved = for molar region)
- convex and blunt working tip
- different widths with conical working end – the different widths are coordinated so that the widest part (black marking) at the working end corresponds to the diameter of the next instrument at the working tip. This ensures continuous expansion of the bone cavity. Moreover, good compression of cancellous bone can be obtained to improve primary stability (especially in the case of soft bone of consistency D3-D4 and D5).

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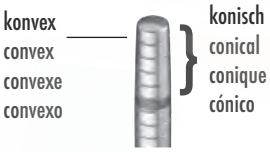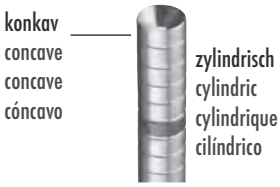
1605-SET


Tray **3029-L-WE**
280 x 190 x 34mm

1605/2.2	1605/2.8	1605/3.5	1605/4.2	1605/4.8	1606/2.2	1606/2.8	1606/3.5	1606/4.2	1606/4.8
Ø2,2mm	Ø2,8mm	Ø3,5mm	Ø4,2mm	Ø4,8mm	Ø2,2mm	Ø2,8mm	Ø3,5mm	Ø4,2mm	Ø4,8mm

OSTEOTOM-SET, Knochenspreizer (Condenser)

OSTEOTOM-SET, Bone spreader + bone condenser / SET D'OSTEOTOMES, Expanseur / SET DE OSTEOTOMOS, Separador de hueso

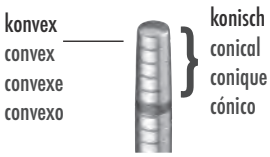1607-SET

Tray **3029-L-WE**
280 x 190 x 34mm

1607/2.2	1607/2.8	1607/3.5	1607/4.2	1607/4.8	1608/2.2	1608/2.8	1608/3.5	1608/4.2	1608/4.8
Ø2,2mm	Ø2,8mm	Ø3,5mm	Ø4,2mm	Ø4,8mm	Ø2,2mm	Ø2,8mm	Ø3,5mm	Ø4,2mm	Ø4,8mm

OSTEOTOM-SET, Knochenmeißel

OSTEOTOM-SET, Bone splitting / SET D'OSTEOTOMES, Ciseaux à os/ SET DE OSTEOTOMOS, Cincel de hueso

doppelseitige Lasermarkierung
double sided laser marking
marquage de laser bilatéral
marcación de láser bilateral

1609-SET

Tray **3029-L-WE**
280 x 190 x 34mm

1609/3
3mm
gerade
straight
droits
rectas

1609/4
4mm
gerade
straight
droits
rectas

1609/6
6mm
gerade
straight
droits
rectas

1609/8
8mm
gerade
straight
droits
rectas


1610/3
3mm
gebogen
curved
courbés
curva

1610/6
6mm
gebogen
curved
courbés
curva

iNNovations

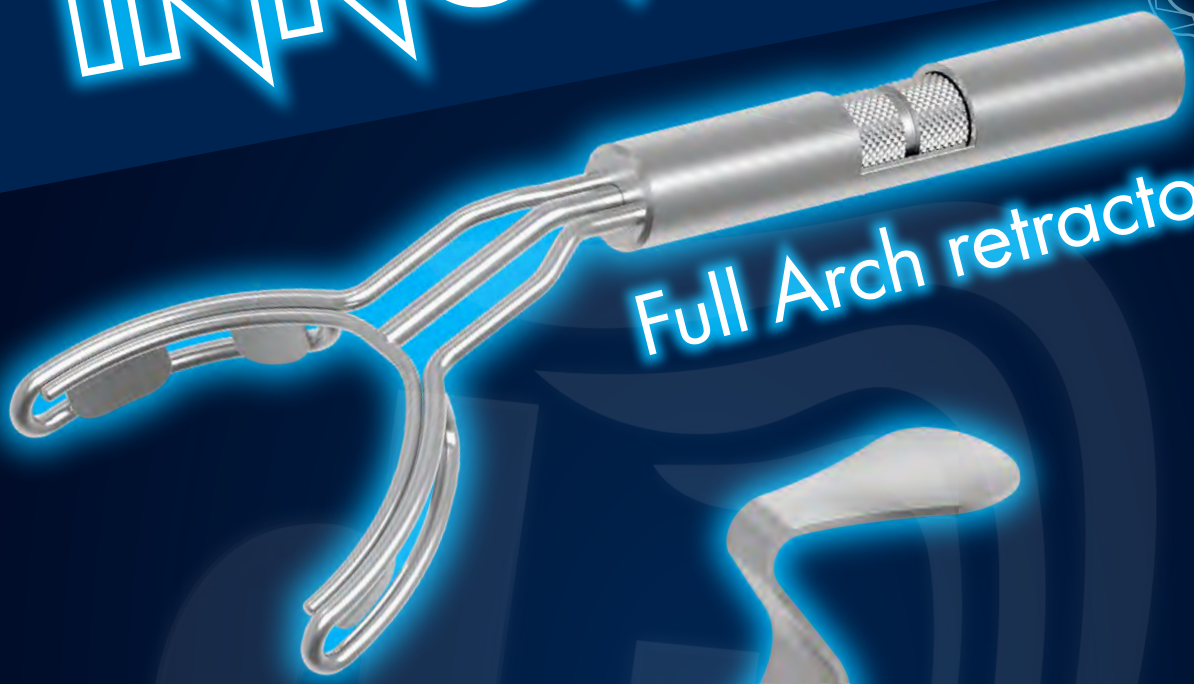Full Arch retractor

Double hook retractor

Zygoma retractor